


Olympus FV3000 ON/OFF Manual

Contact Person:

Lai Ding lding@bwh.harvard.edu

781-901-6116 BTM building 7002M

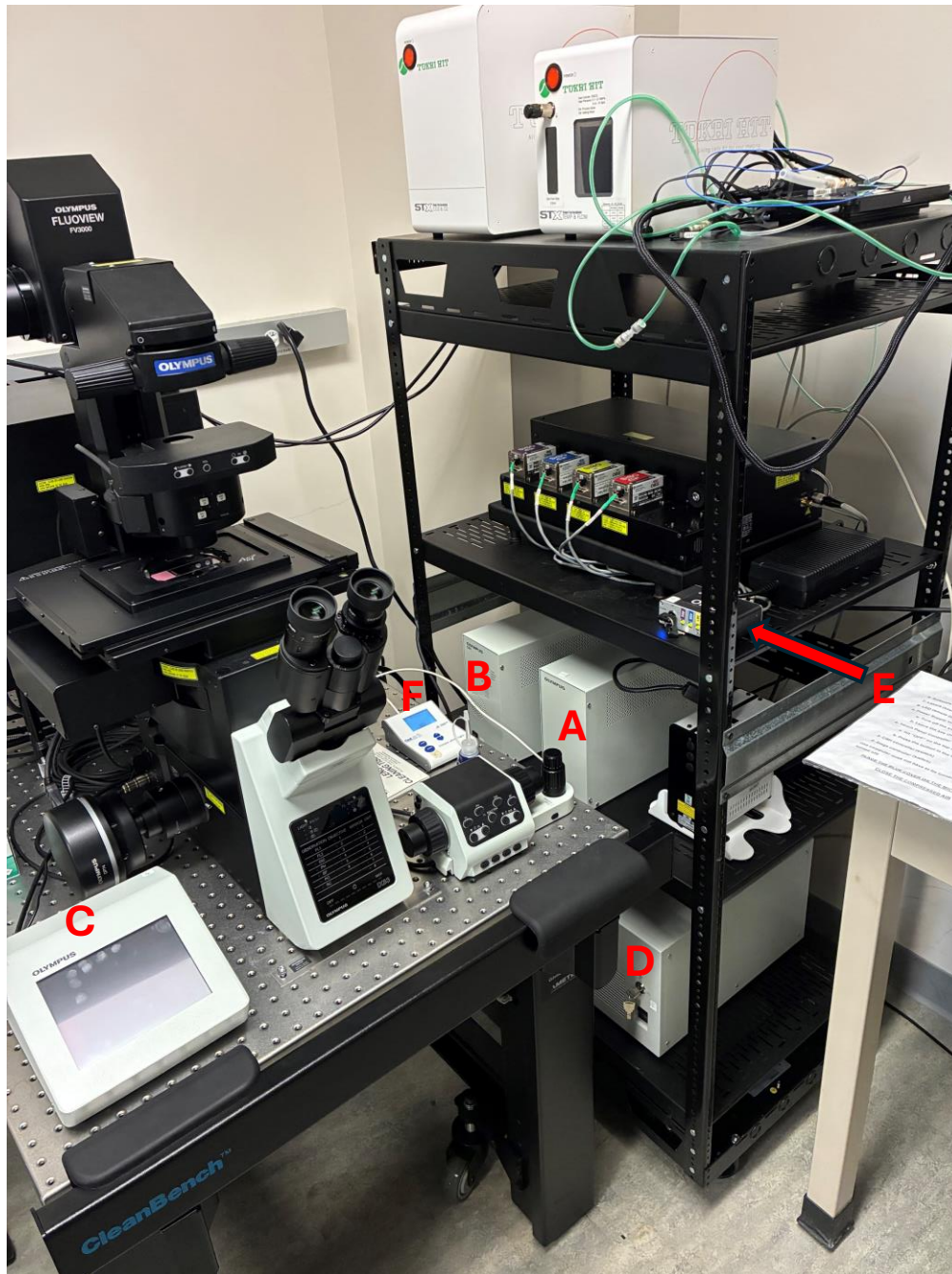

FV3000 system Cold Start up Sequence

A → B → C → D → E → F

A Stage Controller

B CBH Controller

C Touch Panel Controller

D Power Supply Unit (PSU)

E Laser Switch

F LED switch

Simple Switch ON Order (Cold Start)

0. Turn on Computer, Windows has fully started.
- A. Stage Controller Switch ON
- B. CBH Controller Switch ON
Wait for the beep before next step
- C. Touch Panel Controller Switch On (Push button in the back)
Wait for 30 secs until the panel is fully started, then next step
- D. Power Supply Unit Switch ON
Switch ONLY, DO NOT TOUCH THE KEY
- E. Laser key Switch ON.
- F. LED switch ON

Simple Switch OFF Order (If you are the last user)

0. Switch off software FV31S-SW
1. LED Switch OFF **F**
2. Laser key Switch OFF **E**
3. Power Supply Unit Switch OFF **D**
Switch Only, DO NOT touch key
4. Touch Panel Controller Switch OFF
 - a. Hit “OFF” on the top of the screen, then “OK”
 - b. Press the button in the back **C**
5. CBH Controller Switch OFF **B**
6. Stage Controller Switch OFF **A**
7. Computer CAN be switched off or leave on


FV3000 Software ON procedure

1)




Click FluoView icon on desktop

2)

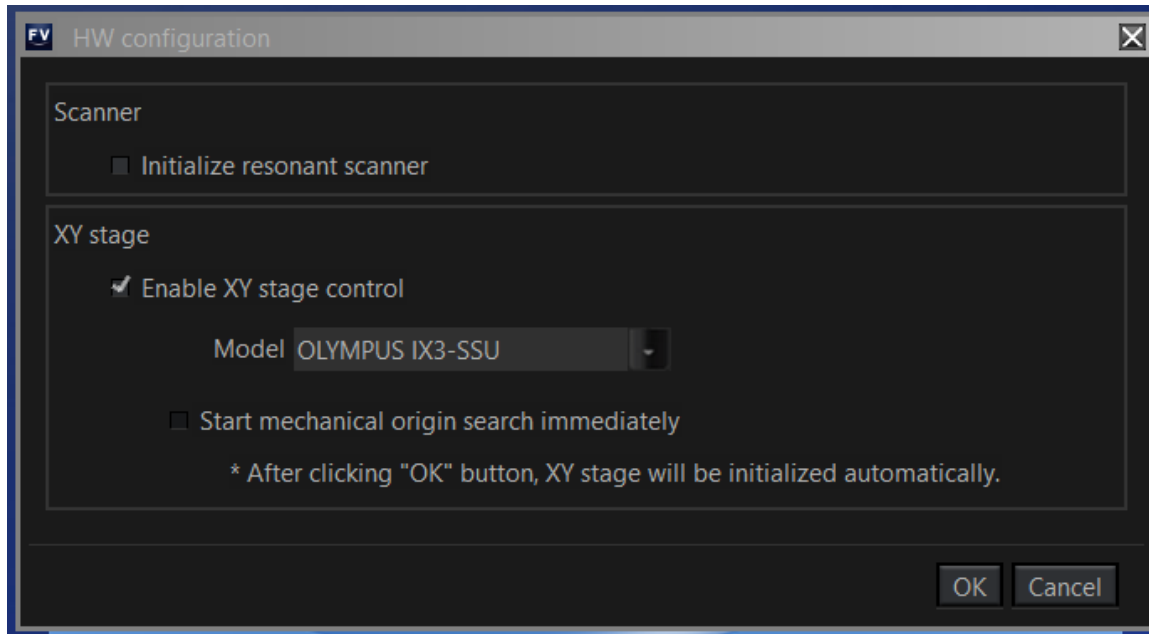

Keep the popup menu as pic on the left, click "OK"

3)

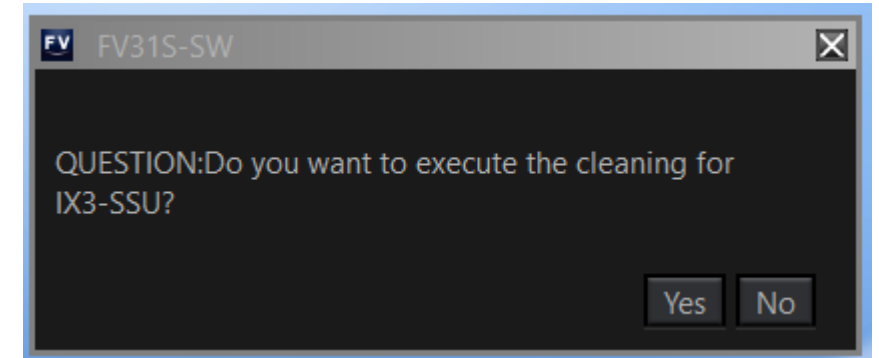

Click NO here !

